

# Flicker Fusion

## iOS User Instructions



3700 Sagamore Parkway North  
Lafayette, IN 47904  
Tel: (765) 423-1505

[info@lafayetteinstrument.com](mailto:info@lafayetteinstrument.com)  
[www.lafayetteinstrument.com](http://www.lafayetteinstrument.com)

## Table of Contents

### Viewing Chamber Hardware

Description	3
Overview	4
Photosensitivity Warning	4
Specifications	5
Parts	5
Operating Instructions	6
Cleaning Instructions	6

### Flicker Fusion App

Pairing Viewing Chamber to App	7
Manage Subjects	8
Add New Subject	8
View and Edit Subjects	8
Task Setup	9
Run Test	9
View Task Results	9
Run Test Interface	10
Automatic Method	12
Adaptive Method	14
Self-Control Method	16
Manual Control Panel	18



### Description

The Lafayette Instrument Flicker Fusion system provides the user with a variety of versatile controls to perform accurate and timely measurements of CFF. Digital circuitry is used to provide extremely accurate frequency generation from 1.0 Hz to 100.0 Hz in 0.1 Hz steps. Unlike previous systems, this generation uses a Bluetooth® connection for maximum convenience. The tablet allows for exact and repeatable frequency settings, and a means for test result scoring and storage. The test options include five modes of operation to cover virtually any test requirement:

- **Simple Flicker:** Stimulus is decreased from steady to the point at which the light changes from steady to flicker.
- **Simple Fusion:** Stimulus is increased from flicker to the point at which the light changes from flicker to steady.
- **Automatic:** Automatically runs test using the Simple Fusion test to determine the fusion frequency and the Simple Flicker test to determine the flicker frequency.
- **Adaptive:** Displays fixed frequencies and allows the test subject to respond to the frequencies as “flicker” or “fusion”. This test zeros in on the CFF and automatically ends when the value is determined.
- **Self-Control:** Allows the test subject to select whether to increase or decrease the frequency. This test is ended by the subject.

The viewing chamber has two lights, one for the left eye and one for the right eye. The light compartments are completely separated, allowing for four stimulus combinations to occur:

- **Left:** Left eye only
- **Right:** Right eye only
- **Coincident:** Left and right simultaneously
- **Alternate:** Left and right alternately

The viewing chamber is constructed to mitigate extraneous factors, which might distort CFF values. The white lamps produce even illumination over a 1” square LED. The stimuli are separated by 2 3/4” center to center with a stimulus to eye distance of 15” and a viewing angle of 1.9 degrees.

### Overview

The Flicker Fusion is designed to test for the critical flicker fusion threshold (CFF). CFF can be used as an index of the temporal resolving power of the human visual system. The threshold of the flicker fusion is determined in the following manner: a light beam is interrupted intermittently by electronic means at a slow rate, causing it to flash or flicker. If the rate of the flicker exceeds a certain rate, the light will appear steady. The rate at which flicker ceases and the light appears steady is the threshold of fusion. The rate of flicker is then decreased from steady to the point at which the light changes from steady to flicker, the threshold flicker. The average of the threshold of fusion and the threshold of flicker is called the Critical Flicker Fusion threshold (CFF).

CFF is sensitive to a number of presentation and observer variables. Presentation variables include stimulus frequency, luminance, and size, color, and contrast, as well as light/dark ratio. Observer variables reported to affect results include body temperature, practice, general physiological state, age, sex, and genetics.

### Photosensitivity Warning

A small percentage of individuals may experience symptoms of photosensitivity when exposed to flashing lights. These symptoms include feelings of nausea, dizziness, migraines, visual distortions, or photosensitivity induced seizures characterized by one or more of the following: light-headedness, altered vision, eye or face twitching, jerking or shaking of arms or legs, disorientation, confusion, loss of consciousness or convulsions that can lead to injury from falling or collision. Symptoms can occur in individuals with no history of seizures or epilepsy. Immediately stop using the Flicker Fusion if you experience any of the above symptoms. Consult a physician before using the Flicker Fusion if you or a relative have a history of seizures or epilepsy.

## Specifications

- **Power Supply**

- AC/DC Wall mount adapter 5V 5W w/ USB 1.1 Cable A Male to B Male 4.92' Shielded
- Computer USB port
- **Frequency:** 1.0 to 100.0 Hz in 0.1 Hz increments with an error of 0.05%
- **Slide Holder:** 2" x 2" 35 mm holder for optional neutral density filters with 0.5% to 50% light transmission
- **Auto Mode Ramp Rates:** Options of 0.5, 1, 2, and 4 Hz per second
- **External Response:** SPDT normally open hand-held switch with 3.5 mm stereo plug
- **Typical Maximum Luminance:** 58 Cd/m<sub>2</sub>
- **Viewing Angle:** 1.9 degrees
- **Light/Dark Ratio:** 1:1
- **Stimulus Color:** White
- **Viewing Chamber Color:** Black
- **Viewing Chamber Mask:** Hobart 770094 Welding Oxy-Acetylene
- **Viewing Chamber:**
  - **Dimensions:** 7" x 19" x 16" (L x W x H)
  - **Weight:** 6.3 lbs

## Parts


- Tablet
- Viewing Chamber
- AC adapter
- Hand-held response switch

### Operating Instructions

1. Connect the black hand-held subject response switch to the external response connector on the rear of the viewing chamber.
2. Plug in wall mount AC adapter via USB port and power on viewing chamber.
3. Enable Bluetooth on tablet and open Flicker Fusion Application
4. Pair viewing chamber with tablet. When the chamber is successfully paired, the wireless connection indicator will light up.

**Note:** For more specific details about the pairing process, see the “Pairing Flicker Fusion viewing chamber to app” section of this manual.

5. Follow app instructions to add users and begin testing.

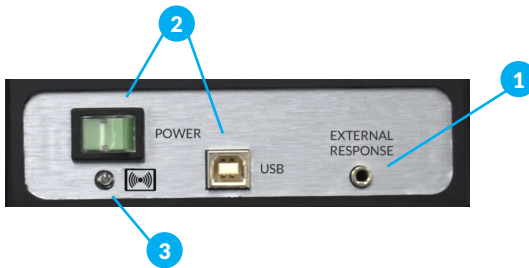

### Cleaning Instructions


When necessary, the mask and viewing chamber can be wiped using a disinfecting cloth or wipe.

### Flicker Fusion App

The Flicker Fusion is controlled by a tablet application. The app allows quick and simple setting up of tests, test scoring and storage of results.

#### Pairing Flicker Fusion viewing chamber to App

The first time the app is loaded or when a new viewing chamber is used, the chamber must be paired to the app. Tap the menu icon [≡] then “Select Device” to view the Select Device interface. Available in this section are a list of currently and previously paired devices as well as a search icon [🔍] to scan for a new viewing chamber.

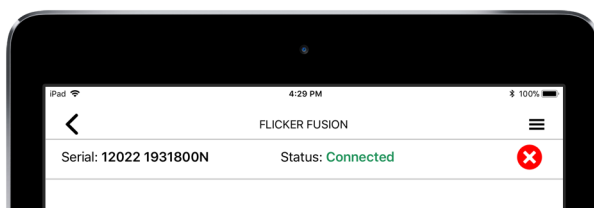




Once a device has been paired, the serial number and status will be displayed below the title bar. A paired device can be disconnected by tapping the [✖] icon.

🔴 **Note:** After disconnecting, the most recently paired device serial number and status will remain in the title area. Quickly reconnect by tapping the [⊕] icon.



### Actions



#### Manage Subjects

Upon opening the application, the user is given the opportunity to manage subjects. From the Manage Subjects interface, tapping the  icon displays the interface for adding a new subject. Tapping on a subject name will display the subject's details. Subject details include previous test and scores, as well as the ability to edit the subject's information and begin a new test. Tapping on the  icon displays the Task Setup interface.

 **Note:** The Manage Subjects interface can also be accessed through the main menu or by tapping the back button.

#### Add New Subject

From the Manage Subjects interface, tapping the  icon displays the interface for adding a new subject. A subject entry includes first and last name and an ID. Tapping on the  icon saves the new subject

 **Note:** Tapping the  icon will save the subject data and continue directly to the Run Test interface.



#### View and Edit Subjects


Subject names are displayed in the Manage Subjects interface. Tapping on an existing subject name displays the Subject Results interface. The Results table is filtered by result date and displays the task title, date, and result. Tapping on an existing test brings up the task result.

 **Note:** From the subject details, tapping on the  icon displays the Edit User interface. Tapping on the  icon displays the run test interface.







### Task Setup



Task setup includes a number of options. These options include task type (Simple flicker, simple fusion, automatic, adaptive, and self-control), stimulus (left, right, coincident, and alternate), luminance (1-100), sweep rate (depends on task type), starting frequency for flicker and fusion (depends on task type and subject), and delay (0-5 seconds). Tapping on the  icon saves the new task setup. Tapping the  icon will run a preview of the chosen task setup. Do not perform testing in preview mode.

 **Note:** For further instructions on testing methods, read the appropriate sections of this manual.


### Run Test

Navigating to the Run Test interface displays the test type, subject name, and various details about how the test will be run. Tapping the  icon begins a new test. Tapping the  icon displays the Task Setup interface.




 **Note:** Tapping the  icon when a test is running will abort the active test. Aborted tests are not scored and cannot be saved.

Once the test has been completed, the examiner can choose to save the test by tapping the  icon or continue without saving by tapping the  icon.


### View Task Results

The View Results interface displays the task title, subject name, date, result, and any chosen test parameters. Tapping on the  icon will delete the chosen test record. Return to the Subject Results interface by tapping the back button.


### Run Test Interface

Regardless of chosen method, all options continue to the Run Test interface. Before the test, only the Run Test  and Test Settings  buttons are enabled while the Abort Test button  is disabled. The example at right is shown during an Automatic test.


1. Display selected test type: Simple Flicker, Simple Fusion, Automatic, Adaptive, or Self Control

 **Note:** For further instructions on testing methods, read the appropriate sections of this manual.

2. Display selected subject (and ID number)
3. Display of current flash rate




 **Note:** Displayed values 4-9 can be edited via the Test Settings interface prior to active testing.

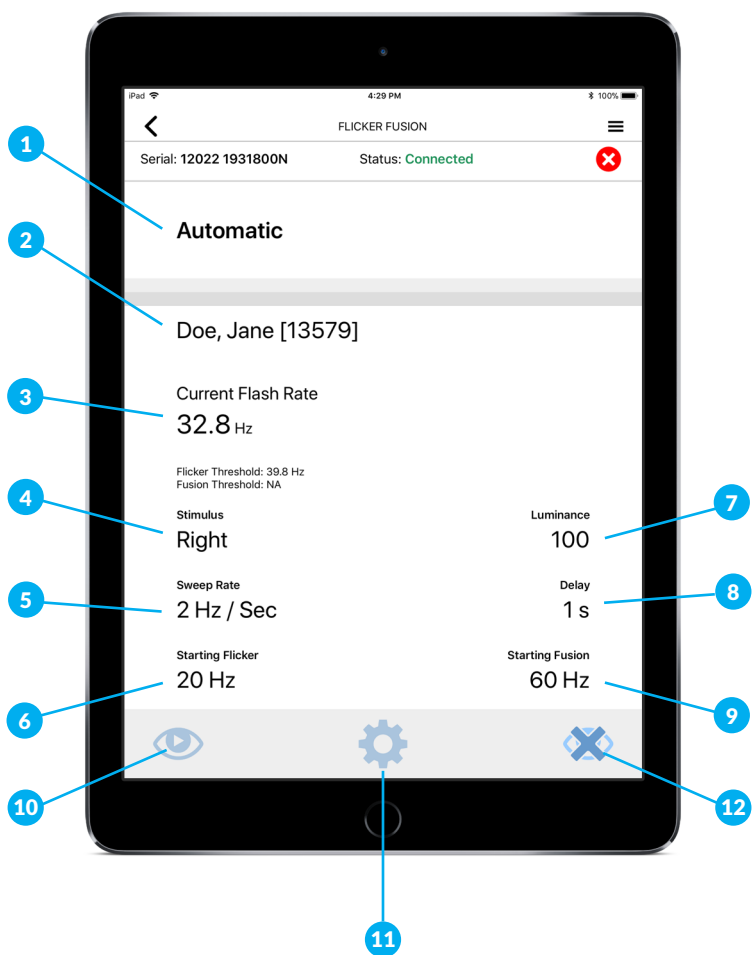
4. Selected Stimulus: Left eye only, right eye only, coincident, or alternate
5. Selected Luminance: Brightness of stimulus
6. Selected Sweep Rate: The rate at which the frequencies will change during the test - shown in Hz/second
7. Selected Delay: Delay between tests
8. Selected Starting Flicker
9. Selected Starting Fusion
10. Run Test button: Begins test - currently disabled due to active test
11. Test Settings button: Displays test settings interface - currently disabled due to active test
12. Abort Test button: Aborts currently active test (enabled due to active test)



### Automatic Method


In Automatic Method, the test is automatically run using the Simple Fusion test to determine the fusion frequency and the Simple Flicker test to determine the flicker frequency. The results of these two tests are averaged to obtain the CFF value. In this method, the subject responds by pressing either side of the hand-held subject response switch.

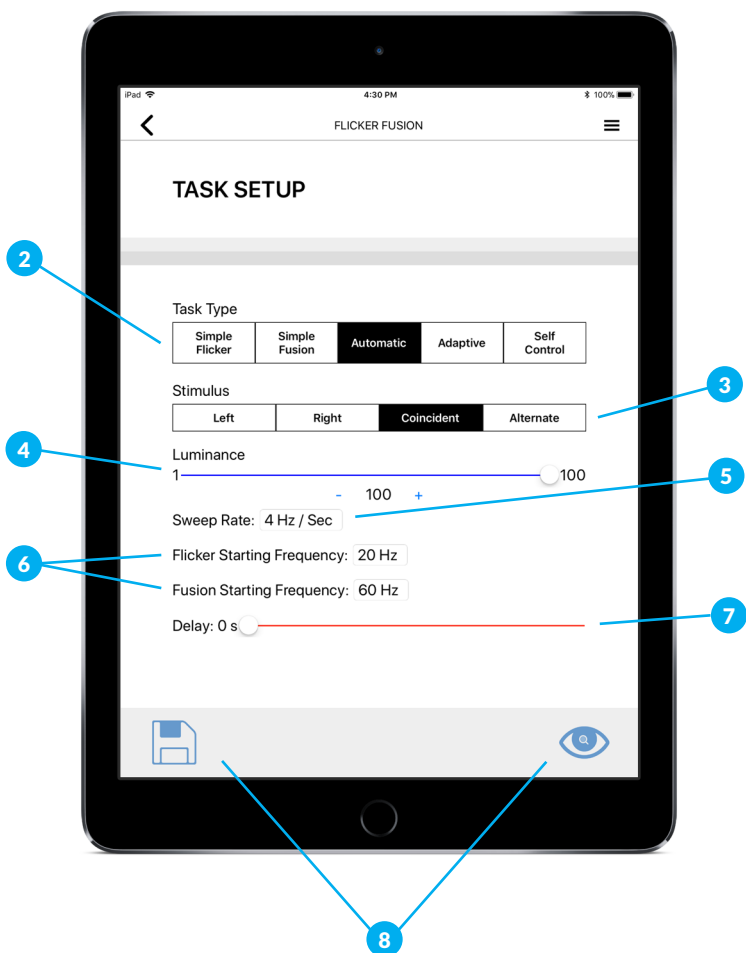


1. From the manage subjects interface, tap the  icon to display the Task Setup interface.
2. Select Automatic Method task from the panel
3. Select the Stimulus (left eye only, right eye only, coincident, or alternate) from the panel. Coincident is selected in the example on the next page.
4. Using the slider, select the desired brightness by adjusting the Luminance (%) setting.
5. The Sweep Rate is the rate at which the frequencies will change during the test in Hz/sec. Select the desired Rate from the appropriate drop down.
6. Select the desired Flicker and Fusion starting frequency by selecting the frequencies from the drop downs. In the example on the next page, the Flicker Starting Frequency is set to 20 Hz and the Fusion Starting Frequency is set to 60 Hz.
7. Using the slider, select the between test Delay Time. Delay Times are measured in seconds within a range of 0 to 5 seconds.
8. Once the appropriate selections have been made, tap the save icon  to save and return to the previous screen OR tap the preview icon  to see a preview of the test using the selected parameters.




### Adaptive Method

In Adaptive Method, the test displays fixed frequencies. Using the hand-held response switch the subject responds to the frequencies as “flicker” or “fusion.” This test zeros in on the CFF and automatically ends when the value is determined.


 **Note:** The half-circle shape with dotted lines extending out on the response switch symbolizes a response of “flicker” while the half-circle shape with solid lines extending out symbolizes a response of “fusion”.

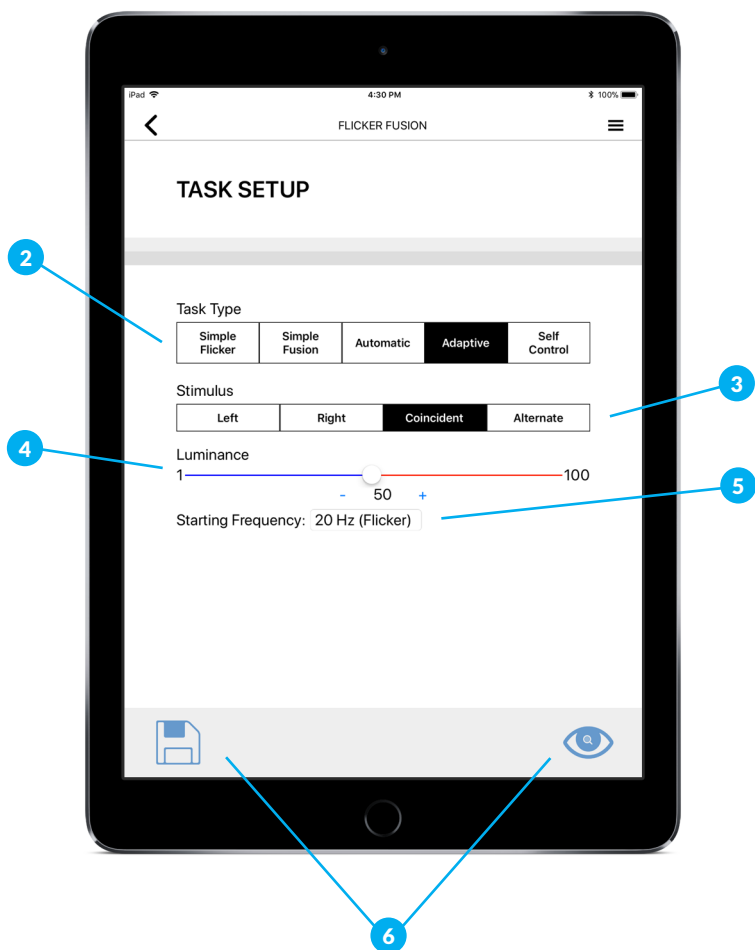


1. From the manage subjects interface, tap the  icon to display the Task Setup interface.
2. Select Adaptive Method task from the panel
3. Select the Stimulus (left eye only, right eye only, coincident, or alternate) from the panel. Coincident is selected in the example on the next page.
4. Using the slider, select the desired brightness by adjusting the Luminance (%) setting.
5. Select the Flicker starting frequency from the drop down. In the example on the next page, the Flicker Starting Frequency is set to 20 Hz.
6. Once the appropriate selections have been made, tap the save icon  to save and return to the previous screen OR tap the preview icon  to see a preview of the test using the selected parameters.



### Self-Control Method

In Self-Control Method, a test allows the test subject to select whether to increase or decrease the frequency. The test is complete when the subject cannot discern whether the frequency is a flicker or a constant. The amount of the increase or decrease becomes smaller as the subject “zeros in” on the CFF. This test is ended by the subject.

 **Note:** The “+” symbolizes a response to increase the frequency while the “-” symbolizes a response to decrease the frequency.

1. From the manage subjects interface, tap the  icon to display the Task Setup interface.
2. Select Self-Control Method task from the panel
3. Select the Stimulus (left eye only, right eye only, coincident, or alternate) from the panel. Right is selected in the example on the next page.
4. Using the slider, select the desired brightness by adjusting the Luminance (%) setting.
5. Using the slider, select the Flicker starting frequency. In the example on the next page, the Flicker Starting Frequency is set to 100.
6. Once the appropriate selections have been made, tap the save icon  to save and return to the previous screen OR tap the preview icon  to see a preview of the test using the selected parameters.






### Manual Control Panel

This interface is used for testing stimulus values by allowing the user to view and select parameters for stimulus, luminance, and flash rate.

1. Tap the menu icon [☰] then “Manual Control Panel” to view the Manual Control Panel interface.
  2. Displays Current Flash Rate
  3. Select the Stimulus (left eye only, right eye only, coincident, or alternate) from the panel. Right is selected in the example on the next page.
  4. Using the slider, select the desired brightness by adjusting the Luminance (%) setting.
  5. Using the slider, select the desired starting Flash Rate setting.
  6. Once the appropriate selections have been made, tap the Flash icon [⚙️] to begin.
- ➡ **Note:** Any changes made on this interface will be immediately reflected by the viewing chamber.
7. Tap the Stop Flash icon [❌] to end.

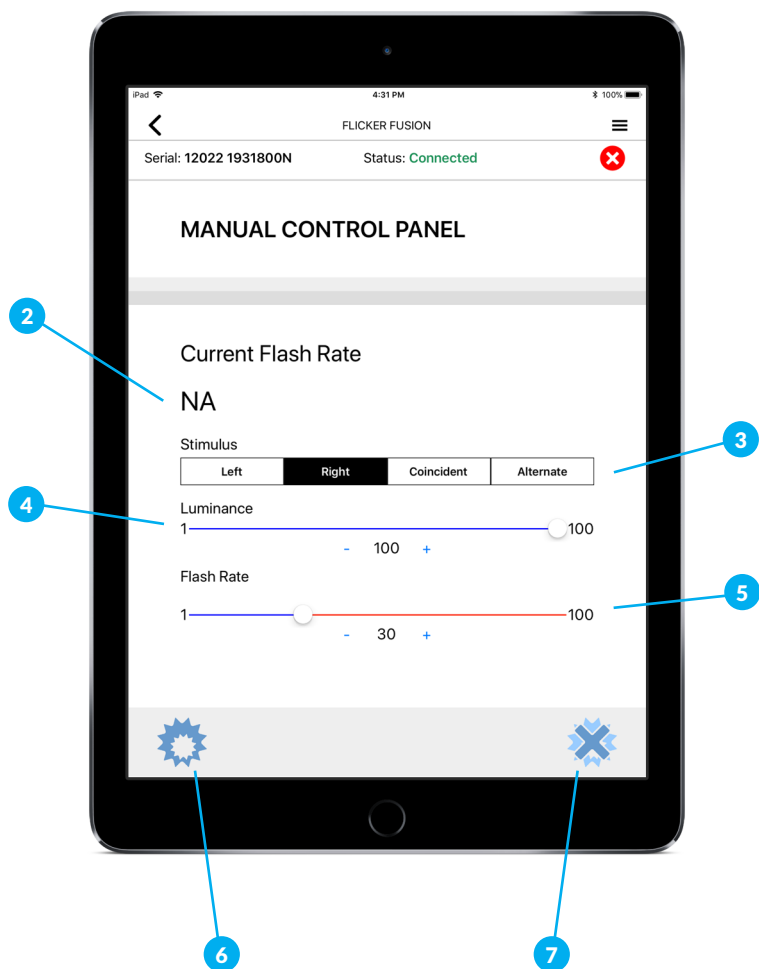

# Terms and Conditions

## LIC Worldwide Headquarters

Toll-Free: (800) 428-7545 (USA only)  
Phone: (765) 423-1505  
Fax: (765) 423-4111  
Email: [sales@lafayetteinstrument.com](mailto:sales@lafayetteinstrument.com)  
[export@lafayetteinstrument.com](mailto:export@lafayetteinstrument.com) (Outside the USA)

## Mailing Address:

Lafayette Instrument Company  
PO Box 5729  
Lafayette, IN 47903, USA

## Lafayette Instrument Europe:

Phone: +44 1509 817700  
Fax: +44 1509 817701  
Email: [eu:sales@lafayetteinstrument.com](mailto:eu:sales@lafayetteinstrument.com)

## Phone, Fax, Email or Mail-in Orders

All orders need to be accompanied by a hard copy of your purchase order. All orders must include the following information:

- Quantity
- Part Number
- Description
- Your purchase order number or method of pre-payment
- Your tax status (include tax-exempt numbers)
- Shipping address for this order
- Billing address for the invoice we'll mail when this order is shipped
- Signature and typed name of person authorized to order these products
- Your telephone number
- Your email address
- Your FAX number

## Domestic Terms

There is a \$50 minimum order. Open accounts can be extended to most recognized businesses. Net amount due 30 days from the date of shipment unless otherwise specified by us. Enclose payment with the order; charge with VISA, MasterCard, American Express, or pay COD. We must have a hard copy of your purchase order by mail, E-mail or fax. Students, individuals and private companies may call for a credit application.

## International Payment Information

There is a \$50 minimum order. Payment must be made in advance by: draft drawn on a major US bank; wire transfers to our account; charge with VISA, MasterCard, American Express, or confirmed irrevocable letter of credit. Proforma invoices will be provided upon request.

## Exports

If ordering instrumentation for use outside the USA, please specify the country of ultimate destination, as well as the power requirements (110V/60Hz or 220V/50Hz). Some model numbers for 220V/50Hz will have a "C" suffix.

## Quotations

Quotations are supplied upon request. Written quotations will include the price of goods, cost of shipping and handling, if requested, and estimated delivery time frame. Quotations are good for 30 days, unless otherwise noted. Following that time, prices are subject to change and will be re-quoted at your request.

## Cancellations

Orders for custom products, custom assemblies or instruments built to customer specifications will be subject to a cancellation penalty of 100%. Payment for up to 100% of the invoice value of custom products may be required in advance. Cancellation for a standard Lafayette Instrument manufactured product once the product has been shipped will normally be assessed a charge of 25% of the invoice value, plus shipping charges. Resell items, like custom products, will be subject to a cancellation penalty of 100%.

## Exchanges and Refunds

Please see the cancellation penalty as described above. No item may be returned without prior authorization of Lafayette Instrument Company and a completed Return Form. A copy of the Return Form or your assigned Return # (you will receive this via email after submitting the form) must be included with the returned goods. The merchandise should be packed well and fully insured. Unopened merchandise may be returned prepaid within thirty (30) days after receipt of the item and in the original shipping carton. Collect shipments will not be accepted. Returned products must be in saleable condition, and credit is subject to inspection of the merchandise.

## Repairs

Instrumentation may not be returned without prior authorization by Lafayette Instrument Company and a completed Return Form. When you complete the

Form, or call Lafayette Instrument, you will receive a Return #: Your Return # number will be good for 30 days. Address the shipment to:

Lafayette Instrument Company  
3700 Saganamore Parkway North  
Lafayette, IN 47904, USA.

Shipments cannot be received at the LIC PO Box. Items should be packed well, insured for full value, and returned along with a copy of the Return Form or the Return #. An estimate of repair will be given prior to completion ONLY if requested in an enclosed cover letter. We must have a completed purchase order by mail or fax, or repair work cannot commence for non-warranty repairs.

## Damaged Goods

Damaged instrumentation should not be returned to Lafayette Instrument prior to a thorough inspection. If a shipment arrives damaged, note damage on delivery bill and have the driver sign it to acknowledge the damage. Contact the delivery service, and they will file an insurance claim. If damage is not detected at the time of delivery, contact the carrier/shipper and request an inspection within 10 days of the original delivery. Please call the Lafayette Instrument Customer Service Department for repair or replacement of the damaged merchandise.

## Limited Warranty

Lafayette Instrument Company warrants equipment manufactured by the company to be free of defects in material and workmanship for a period of one year from the date of shipment, except as provided hereinafter. The original manufacturer's warranty will be honored by Lafayette Instrument for items not manufactured by Lafayette Instrument Company, i.e. resell items. This assumes normal usage under commonly accepted operating parameters and excludes consumable products.

Warranty period for repairs or used instrumentation purchased from Lafayette Instrument is 90 days. Lafayette Instrument Company agrees either to repair or replace, at its sole option and free of part charges to the customer, instrumentation which, under proper and normal conditions of use, proves to be defective within the warranty period. Warranty for any parts of such repaired or replaced instrumentation shall be covered under the same limited warranty and shall have a warranty period of 90 days from the date of shipment or the remainder of the original warranty period whichever is greater. This warranty and remedy are given expressly and in lieu of all other warranties, expressed or implied, of merchantability or fitness for a particular purpose and constitutes the only warranty made by Lafayette Instrument Company.

Lafayette Instrument Company neither assumes nor authorizes any person to assume for it any other liability in connection with the sale, installation, service or use of its instrumentation. Lafayette Instrument Company shall have no liability whatsoever for special, consequential, or punitive damages of any kind from any cause arising out of the sale, installation, service or use of its instrumentation.

All products manufactured by Lafayette Instrument Company are tested and inspected prior to shipment. Upon prompt notification by the Customer, Lafayette Instrument Company will correct any defect in warranted equipment of its manufacture either, at its option, by return of the item to the factory, or shipment of a repaired or replacement part. Lafayette Instrument Company will not be obliged, however, to replace or repair any piece of equipment, which has been abused, improperly installed, altered, damaged, or repaired by others. Defects in equipment do not include decomposition, wear, or damage by chemical action or corrosion, or damage incurred during shipment.

## Limited Obligations Covered by this Warranty

1. In the case of instruments not of Lafayette Instrument Company manufacture, the original manufacturer's warranty applies.
2. Shipping charges under warranty are covered only in one direction. The customer is responsible for shipping charges to the factory if return of the part is required.
3. This warranty does not cover damage to components due to improper installation or use by the customer.
4. Consumable and/or expendable items, including but not limited to electrodes, lights, batteries, fuses, O-rings, gaskets, and tubing, are excluded from warranty.
7. Failure by the customer to perform normal and reasonable maintenance on instruments will void warranty claims.
8. If the original invoice for the instrument is issued to a company that is not the company of the end user, and not an authorized Lafayette Instrument Company distributor, then all requests for warranty must be processed through the company that sold the product to the end user, and not directly to Lafayette Instrument Company.

## Export License

The U.S. Department of Commerce requires an export license for any polygraph system shipment with an ULTIMATE destination other than: Australia, Japan, New Zealand or any NATO Member Countries. It is against U.S. law to ship a Polygraph system to any other country without an export license. If the ultimate destination is not one of the above listed countries, contact us for the required license application forms.